





# Terumo® Spinal Needle



## Patient friendly penetration

- Smooth penetration
- Smooth cutting edge
- Unbendable and unbreakable needle
- Ergonomic, finger-fitting translucent grip-hub
- Easy orientation
- Total colour-coding

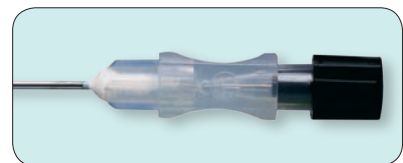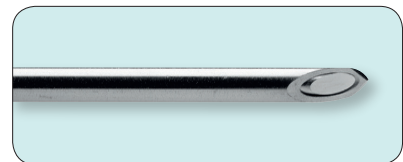

**A smooth penetration with an optimal penetration resistance** thanks to the spinal needle's bevel, which is oriented at a 22° angle. This angle gives you an optimal feel on where exactly the needle point is.

**The smooth and rounded cutting edges** of the cannula and stylet points are perfectly aligned. Due to the precise fit, there is minimal friction between both parts and the stylet can move easily in and out the cannula.

**The needle will not bend or break**, not even in case of an unexpected patient spasm as it is made of special high tensile materials. The danger of an inadvertent detachment is non existing, as both cannula and stylet are firmly affixed to the hub.

**Designed to fit to your fingertips** the grip-hub will ensure that your fingers are not cramped and you can push the needle without straining. The grip-hub will equally allow you to accurately confirm the flow of spinal fluid.

**To easily determine the orientation reading of the needle point**, a TERUMO logo is embossed, aligned with the needle bevel.

**A total colour-coding** is used on the stylet-hub, the blister pack and the shelf package for rapid gauge identification at every stage.



### Terumo® Spinal needle

Product Code	Colour Code		Dimensions of cannula		
			Diameter		Length
SN*1870		Pink	18G	(1.2 mm)	70 mm
SN*1890		Pink	18G	(1.2 mm)	90 mm
SN*1970		Brown	19G	(1.1 mm)	70 mm
SN*1990		Brown	19G	(1.1 mm)	90 mm
SN*2038		Yellow	20G	(0.9 mm)	38 mm
SN*2070		Yellow	20G	(0.9 mm)	70 mm
SN*2090		Yellow	20G	(0.9 mm)	90 mm
SN*2138		Green	21G	(0.8 mm)	38 mm
SN*2170		Green	21G	(0.8 mm)	70 mm
SN*2190		Green	21G	(0.8 mm)	90 mm
SN*2238		Black	22G	(0.7 mm)	38 mm
SN*2270		Black	22G	(0.7 mm)	70 mm
SN*2290		Black	22G	(0.7 mm)	90 mm
SN*2370		Light blue	23G	(0.6 mm)	70 mm
SN*2390		Light blue	23G	(0.6 mm)	90 mm
SN*2590		Blue	25G	(0.5 mm)	90 mm

Sterilized by gamma.  
Unit box: 20 pieces  
Shipping carton: 200 pieces  
Sterilized by ethylene oxide



Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, Terumo Europe provides you with the clinical systems, means and solutions to keep life flowing.



**TERUMO CORPORATION**  
44-1, 2-chome  
Hatagaya  
Shibuya-ku  
Tokyo 151-0072  
JAPAN  
Tel. +81/3 3374.8111  
Fax +81/3 3374.8196

**TERUMO EUROPE N.V.**  
Researchpark Zone 2  
Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.38.12.11  
Fax +32/16.40.02.49

**TERUMO EUROPE N.V. BENELUX SALES DIVISION**  
Researchpark Zone 2  
Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.39.25.80  
Fax +32/16.39.25.99  
THE NETHERLANDS  
Tel. 0800/022.03.96  
Fax 0800/022.04.14

**LABORATOIRES TERUMO FRANCE S.A.**  
1, Parc Ariane  
Bâtiment Uranus  
Rue Hélène Boucher  
78284 Guyancourt Cedex  
FRANCE  
Tel. +33/1 30.96.13.00  
Fax +33/1 30.43.60.85

**TERUMO DEUTSCHLAND GmbH**  
Hauptstrasse 87  
65760 Eschborn  
GERMANY  
Tel. +49/61.96.80.230  
Fax +49/61.96.80.23.200

**TERUMO DEUTSCHLAND GmbH Zweigniederlassung**  
Wettingen  
Bahnhofstrasse 88  
CH-5430 Wettingen  
SWITZERLAND  
Tel. +41/56.419.10.10  
Fax +41/56.419.10.11

**TERUMO EUROPE ESPAÑA SL**  
Avda. Juan Carlos I,  
N° 13 - 7ª Planta  
Edificio Torre La Garena  
28806 Alcalá de Henares  
(Madrid)  
SPAIN  
Tel. general:  
+34/902.10.12.98  
Tel. cust. serv.:  
+34/902.20.31.40  
Fax general:  
+34/902.10.13.58  
Fax cust. serv.:  
+34/902.10.17.57

**TERUMO ITALIA SRL**  
Via Simone Martini  
143/145  
I-00142 Roma  
ITALY  
Tel. +39/06.519.61.420  
Fax +39/06.503.04.07

**TERUMO SWEDEN AB**  
Nya Varvet, Byggnad 90  
Sven Källfelts Gata 18  
SE-426 71 Västra  
Frölunda  
SWEDEN  
Tel. +46/31.748.58.80  
Fax +46/31.748.58.90

**TERUMO DENMARK Filial af TERUMO SWEDEN AB**  
Tel. +45/7020.93.80  
Fax +45/7020.94.80

**TERUMO UK Ltd.**  
Tamesis  
The Causeway  
Egham  
Surrey  
TW20 9AW  
UNITED KINGDOM  
Tel. +44/1784.476.200  
Fax +44/1784.476.201

**TERUMO EUROPE N.V. Moscow Representative office**  
4, Ivana Franko St  
121108 Moscow  
RUSSIA  
Tel. +7/495.234.3303  
Fax. +7/495.234.3303

®: Registered Trademark  
Published by Terumo Europe N.V.



We keep life flowing

Addresses are subject to change, please consult our website.

[www.terumo.com](http://www.terumo.com)  
[www.terumo-europe.com](http://www.terumo-europe.com)